

NATSUME RAT



Product page



Together with the Future of Life Science

“The concept of 3Rs^{*1}” proposed by Russell and Burch in 1959 is currently the international principle on performing animal experiments, and this is an essential principle to follow in the life sciences.

^{*1}Replacement, Reduction, and Refinement

Nazme Co.Ltd. developed and now launch the NATSUME RAT,a silicon mannequin for training animal experimental techniques, taking account of how to address the principle of the 3Rs as a manufacturer of animal experiment devices.

We hope to contribute to a reduction in the number of experimental animals used, and NATSUME RAT is used by those concerned as much as possible as practical efforts for alternatives to animal experiments in medical, pharmaceutical, and veterinary education and the education for sophisticating animal experiment techniques.

【Manufacturer】



NAZME



Our Products

Nazme Co. Ltd.

2-18-6, Yushima, Bunkyo-ku, Tokyo 113-8551 JAPAN

Tel : +81-3-3813-3251 Fax : +81-3-3815-2002

HP : <https://www.nazme.co.jp/>

KN-590

NATSUME RAT

Specification

Length of body	190mm (7½")
Length of tail	180mm (7")
Sex	Male
Body weight	about 200 g
Material	Urethane


Component parts

Main Body of rat model	× 1
Tail	× 2
Instruction manual	× 1

Related Products (optional)

Simulated blood (50 ml)	× 1
Spare tail	× 1

Structure

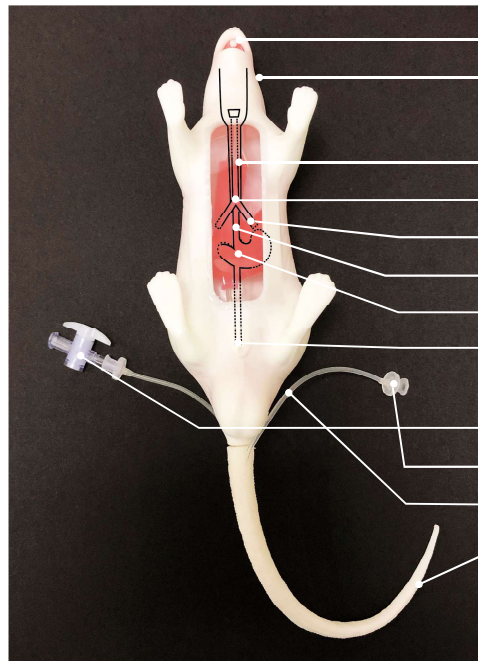

- ① Teeth
- ② Skin
- ③ Trachea
- ④ Bronchial bifurcation
- ⑤ Bronchus
- ⑥ Esophagus
- ⑦ Stomach
- ⑧ Drainage port
- ⑨ Simulated blood cock
- ⑩ Inlet of simulated blood
- ⑪ Tail vein blood vessel tube
- ⑫ Tail


Trainable procedures




Retention /
Holding



Oral administration

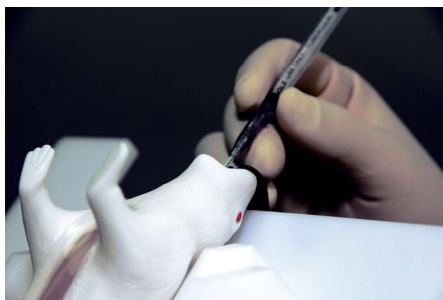

Endotracheal
intubation

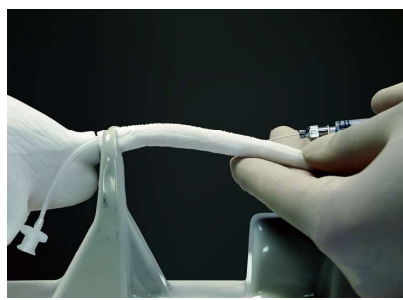

Tail vein injection
Blood collection
through tail vein

- Please, treat NATSUME RAT carefully like a living body. Unnatural handling may cause damage. Refer to the instruction manual for the precautions.
- Please, be aware that the appearance and specifications may change for product improvement without a preliminary announcement.